


# AccuPAP™ fast start guide

With AccuPAP™, you can treat atelectasis effectively and efficiently. Just follow these simple steps.

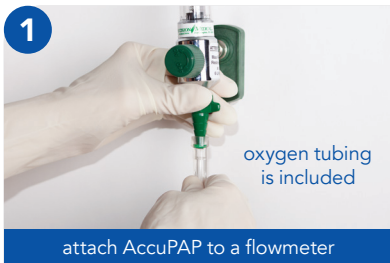

# AccuPAP™ fast start guide

**positive airway pressure therapy has never been easier!**

**adjustable pressure valves** .....  
eliminates the need for a manometer

**patient end with 22mm OD fitting** .....  
can be used with a mask, trach adapter, or during nebulization

**removable mouth piece** .....  
allows for patient's preferred method of treatment delivery

**gas inlet port**  
connects via O<sub>2</sub> tubing (7' included) to wall air or O<sub>2</sub> flow meter

**pressure manometer port**  
for optional pressure monitoring