

# Control-Cric<sup>®</sup>

## Surgical Cricothyrotomy Kit

The Control-Cric is packaged in 3 clear steps, guiding the provider through an otherwise chaotic procedure. Step 1 is the Cric-Knife™, a 10mm dual-sided blade with an integrated sliding tracheal hook. This allows the provider to maintain control of the airway instead of reaching for additional tools and potentially losing the airway. Step 2 is the Cric-Key™, a 5.5mm cuffed tube with an integrated introducer. The introducer gives tactile feedback of the tracheal rings, confirming placement. Step 3 is the neck strap, a flexible silicone strap to slide easily over the tube to maintain position during ventilation.

*Compact 3-Step packaging*

*Rapid deployment to oxygenate/ventilate*

*BVM Wedge included*




### Cric-Knife™

- 10mm dual-sided blade
- Integrated sliding tracheal hook
- Designed to grip like a pencil



### Cric-Key™

- 5.5mm cuffed tube
- Integrated introducer
- Provides tactile feedback of the tracheal rings to avoid false passage

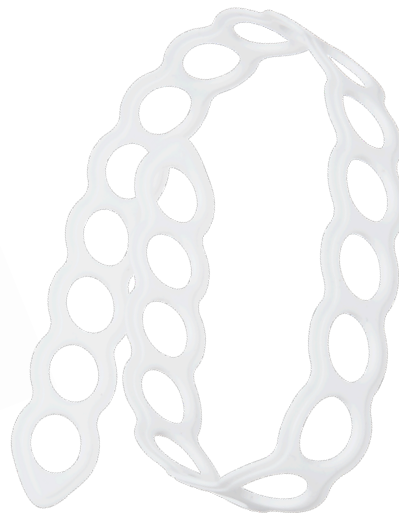

### Neck Strap

- Secure the tube after procedure
- Fits on 22mm connection



### Extension Tube

- Top 15mmF x 22mmM
- Bottom 22mmF
- Swivel Elbow

Part No.	Description	Case/Qty.
313-4567	Control-Cric Kit (Sterile)-Cric-Knife, Cric-Key, Neck Strap, Syringe, BVM Wedge, and Extension Tubing	5
313-4568	Control-Cric Training Kit (Non-Sterile)-Cric-Knife, Cric-Key, Neck Strap, Syringe, BVM Wedge, and Extension Tubing	10

9130F-03-10-2022